

>TruTip[®] Universal Kit


Simplify Your Workflow with a Single Universal Kit for:

- RNA/DNA from Whole Blood and Plasma
- Microbial DNA from Sputum/Sediment
- gDNA from Saliva
- Viral RNA from NPA/NPS

OVERVIEW


Reduce the time of your extraction protocols and simplify your workflow with one Universal Kit for multiple applications. The Akonni Biosystems TruTip Universal Kits offer an innovative solution for efficient isolation of RNA and DNA from challenging clinical sample types for downstream PCR and sequencing applications in just minutes.

- Faster than other methods
- Reduces consumable cost and waste
- Minimal setup for less hands-on time
- High quality nucleic acid

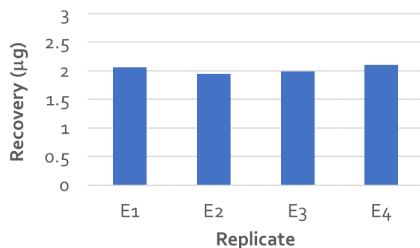


2µg of gDNA recovered from 200µL of whole blood in less than 20 min: 10 min of proK incubation + 9.5 min extraction

Up to 3x Faster Protocol!

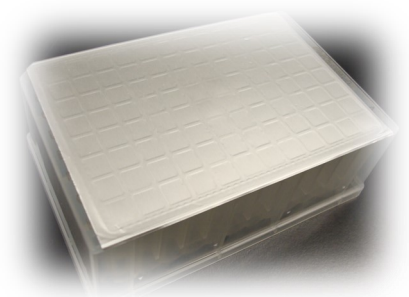


Reduce setup time with all extraction reagents already prepared !


STEP 1: Bind Acid to Matrix

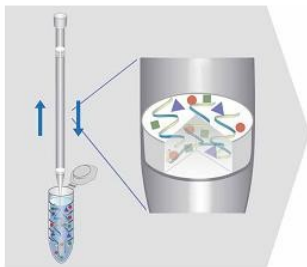


STEP 2: Wash Away Impurities

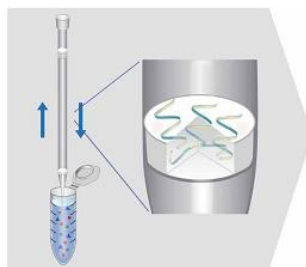

STEP 3: Air Dry

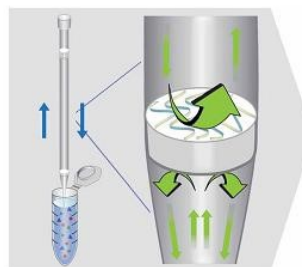

STEP 4: Elute

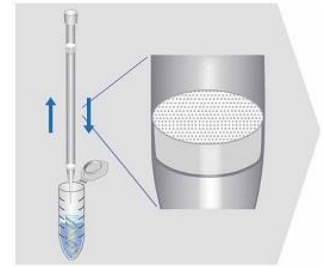

ORDERING INFORMATION

Catalog #
300-30120

Description

TruTip Universal Kit (96 Extractions) For use with the E4 Rainin (Mettler Toledo) Electronic Pipette or Akonni's TruPrepR™ Sample Prep Workstation

For Research Use Only. Not for use in diagnostic procedures. No claim or representation is intended for the diagnosis, prevention, or treatment of a disease. US Patents #7759112, 8,236,501; 8,236,553; 8,574,923; 9,428,746 & additional patents pending. The Akonni logo is a trademark of Akonni Biosystems, Inc. © 2024 Akonni Biosystems. All rights reserved.

Akonni
Biosystems

CORPORATE HEADQUARTERS

400 SAGNER AVE, SUITE 300 FREDERICK, MD 21701 PH: 240-415-9963 FX: 301-205-5785 WWW.AKONNI.COM